

MEDIVATORS™  
**JET PREP™**  
Endoscopic Flushing Device



# MEDIVATORS™ JET PREP™

Endoscopic Flushing Device




High powered irrigation and suction to clear and remove various GI obstructions and increase visualization

Studies have shown 27% of colonoscopies have suboptimal bowel preparation<sup>1</sup>, with the right colon being especially problematic. This subsequently increases the potential for missed adenomas. With the JET PREP Endoscopic Flushing Device you can increase your visualization during diagnostic and therapeutic procedures.

**BEFORE**






**JET PREP IRRIGATION**

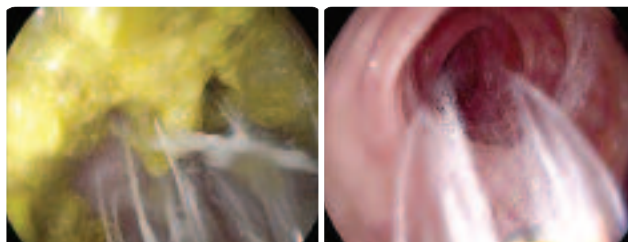

**AFTER**



**SERRATED LESION DETECTION**

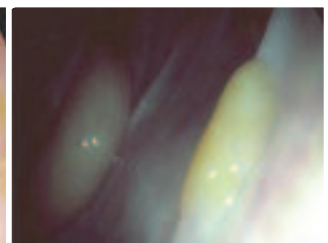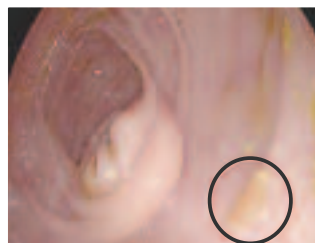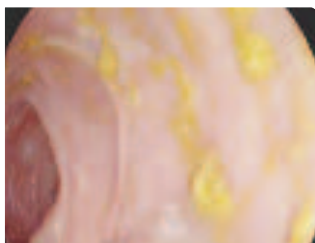

**BEFORE**

**AFTER**

1. Cohen L., Tennyson C., Bowel Preparation for Colonoscopy: Maximizing Efficacy, Minimizing Risk. Gastroenterology & Endoscopy News. 2011

# Tracking ADR? You Can't Find What You Can't See

The JET PREP™ Endoscopic Flushing Device provides an immediate solution for sub-optimal bowel prep, GI bleeding and food impaction.

## JET PREP Irrigation



- High-powered spray improves visualization and breaks down food impactions

## Suction



- Uninterrupted suction capability
- Helps identify source of bleeding through simultaneous flushing and evacuation

## Channel Clearing



- Reduces clogging of working channel via backflush feature

## Powerful Irrigation Solution

### ENDO STRATUS® Irrigation Pump

- High performance flow
- Heated water option

### ENDOGATOR® Disposable Irrigation Tubing

- 24 hour use without the need to reprocess
- Improve efficiencies and infection prevention practices



## JET PREP™ Endoscopic Flushing Device Scope Connections



OLYMPUS® Clip Channel Connector



PENTAX® Luer Channel Connector



FUJINON® Clip Channel Connector

### JET PREP Flushing Device SKU Compatibility

SKU	Model Type	Scope Type	Channel Diameter (mm)	Scope Length (mm)
JP30022	OLYMPUS	Colonoscope - Long	3.7	1680
JP30016	OLYMPUS	Colonoscope - Short	3.7	1330
JP30069	OLYMPUS	Pediatric Colonoscope - Long	3.2	1680
JP30068	OLYMPUS	Pediatric Colonoscope - Short	3.2	1330
JP30050	OLYMPUS	Gastroscope	3.7	1030
JP30064	PENTAX	Colonoscope - Long	4.2	1700
JP30063	PENTAX	Colonoscope - Short	4.2	1500
JP30024	PENTAX	Pediatric Colonoscope - Long	3.8	1700
JP30023	PENTAX	Pediatric Colonoscope - Short	3.8	1500
JP30051	PENTAX	Gastroscope	3.8	1050
JP30052	FUJINON	Colonoscope - Long	3.8	1690

JET PREP™ is a trademark of Jet Prep Ltd. MEDIVATORS™ is a trademark of MEDIVATORS Inc. ENDOGATOR® and ENDO STRATUS® are registered trademarks of MEDIVATORS Inc. OLYMPUS® is a registered trademark of Olympus Corporation. PENTAX® is a registered trademark of Hoya Corporation. FUJINON® is a registered trademark of FUJIFILM Corporation.



MEDIVATORS Inc. Headquarters  
14605 28th Avenue North  
Minneapolis, MN  
55447-4822 USA  
Tel: +1.800.444.4729  
Fax: +1.763.553.3387

MEDIVATORS BV  
Sourethweg 11  
6422 PC Heerlen  
The Netherlands  
Tel: +31.45.5.471.471  
Fax: +31.45.5.429.695

MEDIVATORS Asia/Pacific Pte. Ltd.  
1 A International Business Park  
#05-01 Singapore 609933  
Tel: +65.6227.9698  
Fax: +65.6225.6848

MEDIVATORS Beijing Representative Office  
Room 708, 7th Floor Kaiheng Center, Block B  
No. 2 Chaoyangmennei Street  
Dongcheng District, Beijing China 100010  
Tel: +8610.6567.8446  
Fax: +8610.6567.8445

[www.medivators.com](http://www.medivators.com)

© 2014 MEDIVATORS Inc. P/N 50098-472A